

KEEP THIS USER GUIDE FOR FUTURE REFERENCE

Always retain your proof of purchase in case of warranty service.


Keep all packaging materials out of reach of children to avoid risk of suffocation.

Oricom MNA01 Nasal Aspirator is included as a medical device on the Australian Register of Therapeutic Goods (ARTG No. 456552) and is intended to enable clearing of excessive mucus from nasal passages of an infant/child to facilitate easier breathing. Intended for domestic use only. ALWAYS READ THE USER GUIDE AND USE ONLY AS DIRECTED. IF SYMPTOMS PERSIST SEE YOUR DOCTOR/HEALTHCARE PROFESSIONAL. Use this product only for the intended purpose as described in this user guide.

CAUTION!

- This product is designed to support your child when they have a cold and is not a medicinal product. If your child is sick, consult a doctor.
- Never blow into the mouthpiece. Not following the instructions can endanger your child's health!
- If nasal secretion is too tough or dry, try to soften the blockage with an isotonic saline solution, spray or drops. In case of doubt or uncertainty, please contact your pediatrician.

DELIVERABLES & PART DESCRIPTION

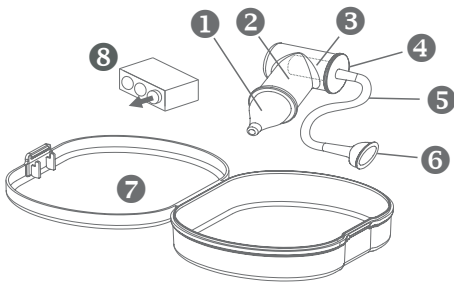

1. Soft Tip
2. Collecting Container
3. Hygiene Filter (Can be washed out)
4. Hose Connector
5. Hose
6. Mouthpiece
7. Storage Box
8. 3pcs Hygiene Filter

SAFETY

- Do not use if there are any signs of damage.
- Do not use the Nasal Aspirator for nose bleeds or nostril wounds.
Do not use on the ears.
- Only for use by adults.
- This product is not a toy, keep out of the reach of children.
- Make sure all parts are connected securely to one another before use.

INSTRUCTIONS FOR USE

1. Insert a new or cleaned filter (3) into the hose connection(4).
2. Connect the hose connection (4) to the collection container (2) and the silicone tip (1).
3. Place the hose (5) onto the hose mount of the collection container.
4. Place the mouthpiece (6) into the end of the hose (5).
5. Place the soft silicone tip tightly against your child's nasal opening.
6. Draw on the mouthpiece until the nasal discharge ends up in the collection container.
7. You can remove the hose connection and replace the filter if desired. The filter prevents the transfer of bacteria and is washable.

 **Place the Nasal Aspirator close to the opening of the nostril. Do not push it into the nostril!**

CLEANING

1. Disassemble the nasal aspirator directly into individual parts after each use.
2. Wash all parts (including the filter) immediately and carefully with warm water and mild detergent.
3. Rinse the clean parts with fresh water.
4. Let the parts dry completely.
5. Clean only in the top basket of the dishwasher or as described in points 2~4.
6. When not in use for long periods, we recommend that you clean the nasal aspirator again and sterilise it before the next use.

 **CAUTION:** DO NOT place this product in a microwave.

STORAGE

- We recommend that you store the product out of direct sunlight.

Express Warranty (Australia)

This Express Warranty is provided by Oricom International Pty Ltd ABN 46 086 116 369, Unit 1, 4 Sovereign Place, South Windsor NSW 2756, herein after referred to as "Oricom".

Oricom warrants that the product is free from defects in materials or workmanship during the Express Warranty Period. This Express Warranty does not extend to any product from which the serial number has been removed or was purchased outside of Australia.

The benefits of this Express Warranty are in addition to other rights and remedies you may have under the Australian Consumer Law. Our goods come with guarantees that cannot be excluded under the Australian Consumer Law. You are entitled to a replacement or refund for a major failure and compensation for any other reasonably foreseeable loss or damage. You are also entitled to have the goods repaired or replaced if the goods fail to be of acceptable quality and the failure does not amount to a major failure. In the event of a minor failure, Oricom reserves the right to choose to repair or replace the product.

The Express Warranty Period will be a period of 12 months beginning on the date of purchase of the product evidenced by your dated sales receipt. You are required to provide proof of purchase as a condition of receiving Express Warranty services.

You are entitled to a replacement product or repair of the product at our discretion according to the terms and conditions of this document if your product is found to be faulty within the Express Warranty Period. This Express Warranty extends to the original purchaser only and is not transferable.

Products distributed by Oricom are manufactured using new materials or new and used materials equivalent to new in performance and reliability. Spare parts may be new or equivalent to new. Spare parts are warranted to be free from defects in material or workmanship for thirty (30) days or for the remainder of the Express Warranty Period of the Oricom branded product in which they are installed, whichever is longer. During the Express Warranty Period, Oricom will where possible repair and if not replace the faulty product or part thereof. All component parts removed under this Express Warranty become the property of Oricom. In the unlikely event that your Oricom product has a recurring failure, Oricom may always, subject to the Competition and Consumer Act 2010, at its discretion, elect to provide you with a replacement product of its choosing that is at least equivalent to your product in performance.

No change to the conditions of this Express Warranty is valid unless it is made in writing and signed by an authorised representative of Oricom.

Oricom will not be liable under this Express Warranty, and to the extent permitted by law will not be liable for any defect, loss, damage or injury arising out of or in connection with a:

1. Failure by you to adhere to the warnings and follow the instructions set out in this user guide for the proper installation and use of the product;
2. Willful misconduct or deliberate misuse by you of the product;
3. Any external cause beyond our control, including but not limited to power failure, lightning or over voltage; or
4. Modification to the product or services carried out on the product by anyone other than Oricom or Oricom's authorised service provider.

How to make a claim under your Express Warranty in Australia

Oricom has a simple warranty process for you to follow:

- Please call or email our Customer Support Team, 02 4574 8888 or support@oricom.com.au.
- A Customer Support Team member will verify after troubleshooting with you if your product qualifies under warranty. If so, they will give you a Product Return Authorisation number.
- We will then email a Return Authorisation form and a Repair Notice (if necessary), together with instructions on how to return the product for warranty service.

Please note that if a Customer Support Team member advises that your product does not qualify for return, this warranty does not apply to your product. Products that are authorised to be returned to Oricom in Australia must include all of the following:

- A completed Return Authorisation form
- A copy of your Proof of Purchase (please keep your original copy)
- The faulty product, including all accessories.

Send the approved returns to:

Oricom International Pty Ltd
Locked Bag 658
South Windsor NSW 2756 Australia

Please note that this Express Warranty excludes expenses incurred by you in returning any faulty product to us. You must arrange and pay any expenses incurred (including postage, delivery, freight, transportation or insurance of the product) to return the faulty product to us, however, we will arrange delivery of the repaired or replaced faulty product to you.

Important Information - Repair Notice

Please be aware that the repair of your products may result in the loss of any user-generated data (such as stored telephone numbers, text messages and contact information). Please ensure that you have made a copy of any data saved on your product before sending for repair. Please also be aware that products presented for repair may be replaced by refurbished products or parts of the same type rather than being repaired.

Contact Details for Oricom support and Warranty claims in Australia

Oricom International Pty Ltd
Unit 1, 4 Sovereign Place
South Windsor, NSW 2756
Australia

Email: support@oricom.com.au
Phone: (02) 4574 8888
(Monday to Friday 8am to 6pm AEST)
Web: www.oricom.com.au